ww[®]
White Way

WW-Sonic

SONIC TOOTHBRUSH :: SONIC-ZAHNBÜRSTE :: SONICZNA SZCZOTECZKA DO ZĘBÓW :: SONICKÝ ZUBNÍ KARTÁČEK

EN DE PL CZ

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USER MANUAL

Sonic Toothbrush

Thank you for buying WhiteWay product.

The sonic technology of WW-Sonic toothbrush provides optimum cleaning of the tooth surfaces and interdental areas inaccessible to traditional toothbrushes. It is equipped with a sound wave generator that generates 31,000 to 48,000 sonic vibrations per minute. It enables the bristles effectively remove plaque, tartar, stains, but also gum massage returning your teeth white and healthy.

Please read these instructions before using.

DANGERS


- Do not immerse totally in water or in any liquids.
- Do not use the product while bathing.
- Do not place or store the product where it can fall into tub or sink. Do not reach for the product that has fallen into water. Unplug immediately.
- Do not disassemble, repair or remodel the main unit. Warranty is void under unauthorized attempt to open the main unit.

WARNINGS

- Toddlers or people who are not able to operate should be supervised by adults.
- Use this product only for its intended use as described in this manual.
- Do not use attachments not recommended by the manufacturer.
- Do not operate any appliance with a damaged cord or plug. Contact authorized service facility for examination.
- Do not plug or unplug the power cord into the electrical outlet that has a voltage different than that specified by the producer.
- Do not use the charging base near heated surfaces or outdoors.
- Do not immerse the charging base in water.
- Do not expose the product to high fluctuation of temperature, direct effect of the sun's rays; do not drop and avoid sudden shocks.
- Never insert any object into any opening.
- Do not use damaged brush heads.
- Never use while sleeping or drowsy.
- Do not put any metal articles on the charging base.
- Oral hemorrhage may occur when use the toothbrush at the first time. If the oral hemorrhage lasts for more than a week, please consult the doctor for suggestion.
- Please consult the doctor for suggestion before use the toothbrush if your body is equipped with heart pacemaker or other electronic devices.

Keep this instruction manual !

CONFIGURATION

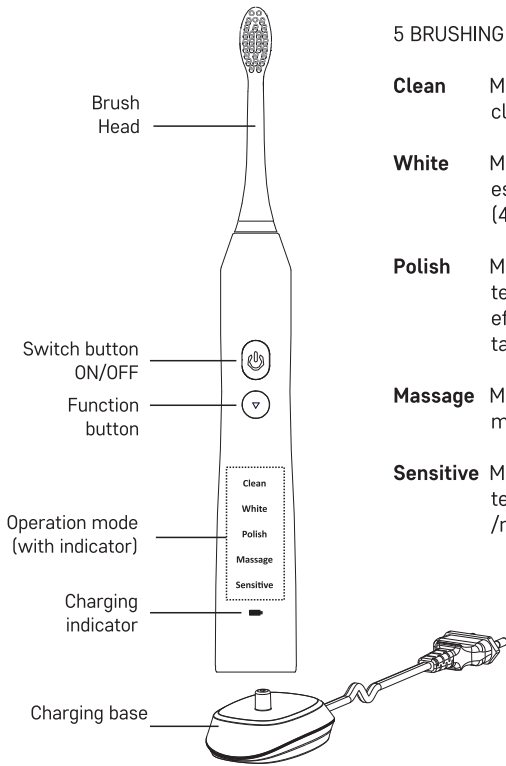
5 BRUSHING MODES:

- Clean** Mode for daily teeth and gums cleaning (31 000 strokes /min.)
- White** Mode for removing surface stains, especially of coffee and cigarettes. (41 000 strokes / min.)
- Polish** Mode for brightening and polishing teeth. With the highest strokes effectively removes stains and tartar. (48 000 strokes / min.)
- Massage** Mode for cleaning and gum massage (31 000 strokes /min.)
- Sensitive** Mode for gentle care of sensitive teeth and gums (31 000 strokes /min.)


CHARGE FOR ABOUT 12 HOURS BEFORE FIRST USE.

1. Place the charging base on a platform, then plug into power supply.
2. Put the handle on the charging base. When charging, the charging indicator illuminates in red while all the mode indicators flash in blue in sequence. The charging indicator turns to blue and all the mode indicators turn off when the sonic toothbrush is fully charged.


OPERATING YOUR SONIC TOOTHBRUSH

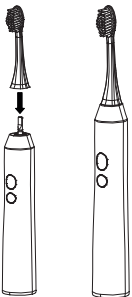
1. Install the brush head onto the handle.

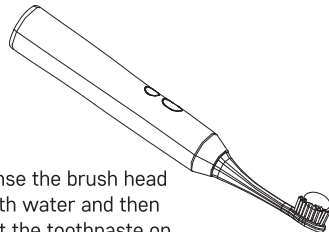2. Rinse the brush head with water and then put the toothpaste on the brush head.

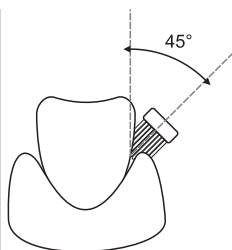3. Put the brush head into mouth.

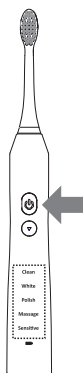4. Press down the switch button, and then press the function button to select the operation mode based on your demand. There are 5 operation modes including "Clean", "White", "Polish", "Massage" and "Sensitive" for your selection. After selecting, its corresponding indicator will illuminate, then the toothbrush starts to work.

5. Gently move the brush on each tooth surface. Do not press the toothbrush against your teeth too firmly.
6. The toothbrush will stop vibrating once about every 30 seconds to remind you to brush other sections of your mouth, and it will shut off automatically after running for 2 minutes. The corresponding indicator of selected operation mode will extinguish after finish working.

NOTE: Press the ON/OFF button once if need to make the toothbrush stop working during brushing process.

7. After cleaning teeth, put the brush head under the running water. If the toothbrush runs out of power, put it back to the charging base.

LOW-VOLTAGE PROTECTION

When the voltage of sonic toothbrush is too low, the charging indicator will flash in red and the toothbrush cannot be operated so as to protect the battery being discharged excessively. At that time, you need to recharge the toothbrush.

CLEANING AND MAINTENANCE

Disconnect the power supply before cleaning. Never immerse the charging base into water. Rinse the handle and the brush head. Use only a soft and dry cloth to clean the charger.

NOTE: Never use a poisonous cleaner to clean all parts of the toothbrush. Do not clean brush head, charging base or handle in the dishwasher.

NOTE: Replace the brush head every 3 months or sooner if signs of wear appear.

WASTE DISPOSAL

Used up product shall be supplied to a waste collection point. Contains ingredients hazardous to the environment. Correct disposal of the device helps secure precious resources and avoid negative impact on health and the environment, which can be endangered by incorrect waste handling. If you have any doubts where to dispose of your device, contact Diagnosis S.A. or your local distributor.

USED SYMBOLS:

Manufacturer

Direct current

Class II insulation protection

Read the operating instructions before use

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EN

Sonic strokes
Power

31000-48000 per minute
Li-ion 800mAh inductive charging provides 180 minutes working time
worldwide voltage 100-240V
Working power 1.4W, standby power 0,5W
≤ 65dB
50x75x247mm
IPX7
2 min.
30 sec.
18 - 25°C; 10-95% RH
-20 - 60°C; 10-95% RH
2 years
sonic toothbrush, two brush heads, inductive charger, user manual

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WARRANTY TERMS

EN

- Diagnosis S.A. grants a warranty: 24 months for Toothbrush WW-SONIC (excluding brush heads).
 - Hardware defects revealed during the warranty period will be rectified free of charge within 21 days.
 - The term runs from the date of delivery of the equipment to the service center.
- The purchaser shall be entitled to replace the equipment for a new one, free of defects, when:
 - The repair has not been made within the time limit set in item 1
 - An authorized service center found an irreparable manufacturing defect
 - During the warranty period, 4 repairs were effected, and the equipment still shows defects that prevent its use in accordance with its intended purpose.
 - The concept of repair shall not include operations related to equipment check and cleaning.
- The warranty shall not cover: batteries, products with illegible or damaged serial number, damage due to the operation and storage inconsistent with the user manual, ingress of liquids or foreign bodies, overvoltage of mains, repairs by unauthorized persons and random events.
- Faulty equipment should be delivered by the buyer to the address of the main service center.
- The warranty for the sold consumer goods shall not exclude, restrict or suspend the powers of the buyer resulting from non-conformity of the goods with the contract.
- The only basis for the warranty rights shall be the warranty card with the date of sale, stamp and signature of the salesperson. If the card is not completed, filled in wrongly, with traces of corrections and entries made by unauthorized persons, illegible as a result of damage - it shall be invalid.

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The buyer should provide the equipment in a properly protected package. Damages caused by insufficiency or unsuitability of packing are excluded from the warranty. Diagnosis S.A. is also not liable for defects that resulted from the damaging caused by post or courier companies.

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- Diagnosis S.A. udziela gwarancji: 24 miesiące na szczoteczkę WW-SONIC (z wyłączeniem końcówek)
- Diagnosis S.A. poskytuje garanci 24 měsíců na sonický zubní kartáček WW-SONIC (bez koncovek)